


MetPrep Preparation Procedure – No 195

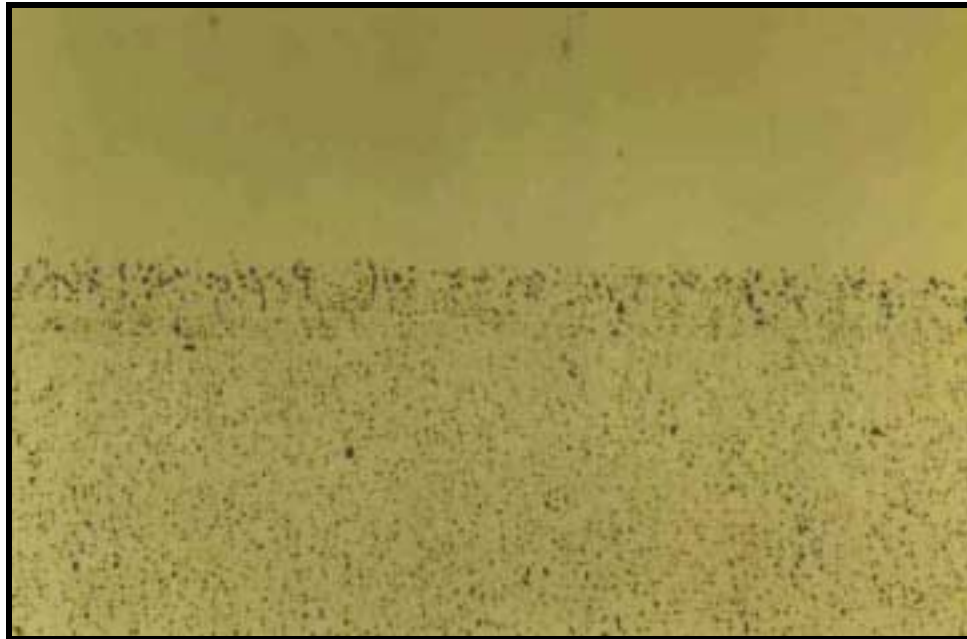

Silver – 10% Cadmium Wire

	Surface	Abrasive	Pressure		Speed – Dir	Time
Primary Grinding Stage	Paper	P240g SiC	Psi	N	150 – Comp	Until Planar
			5	25		

	Surface	Abrasive	Pressure		Speed – Dir	Time
Additional Grinding Stages	Planocloth H	9 µm (OB) PCD Diamond	Psi	N	250 – Comp	4 mins
			5	25		
	Durasilk	3 µm (OB) PCD Diamond	5	25	250 – Comp	4 mins
	Durasilk	1 µm (WB) PCD Diamond & 0.06µm Silco	5	25	100 – Comp	3 mins

	Surface	Abrasive	Pressure		Speed – Dir	Time
Polishing Stage	Multicloth	0.06µm Silco	Psi	N	80 – Comp	2 mins
			3	15		

MetPrep Ltd, Curriers Close, Charter Avenue, Coventry. CV4 8AW
Tel: +44 (0) 24 7642 1222 Fax: +44 (0) 24 7642 1192
Website: www.metprep.co.uk Email: sales@metprep.co.uk